



WiScope® Image System

◎ CYSTOSCOPE SPECIFICATIONS

Field of View	100°
Direction of View	Forward
Depth of Field	2 – 50mm
Distal Tip Diameter	9.0Fr
Insertion Tube Outer Diameter	13.8Fr
Working Channel Diameter	6.6Fr
Angulation Range	180° Up / 180° Down
Working Length of Shaft	380mm
Total Length	615mm

◎ IMAGE SYSTEM SPECIFICATIONS

Voltage	AC 100-240V
Frequency	50/60Hz
Key Functions	Brightness Increase
	Brightness Decrease
	Zoom
	Image Freezing
Video Output	CVBS/HDMI Output to Monitor
	USB Output to Computer
Dimensions	Width 32.5cm × Height 6.5cm × Depth 30.5cm
Weight	2.7kg

CONTACT

OTU Medical Inc.
 1828 Bering Dr, San Jose, CA 95112, USA
 www.otumed.com
 sales@otumed.com



WiScope® Single-Use
 Digital Flexible Cystoscope




- No cross-infection
- No sterilization
- No maintenance

Bring smart views to One Time Use (OTU) surgical devices

MODELS

OTU-C380SR: Right-Hand Standard Deflection
OTU-C380RR: Right-Hand Reverse Deflection
OTU-C380SL: Left-Hand Standard Deflection
OTU-C380RL: Left-Hand Reverse Deflection

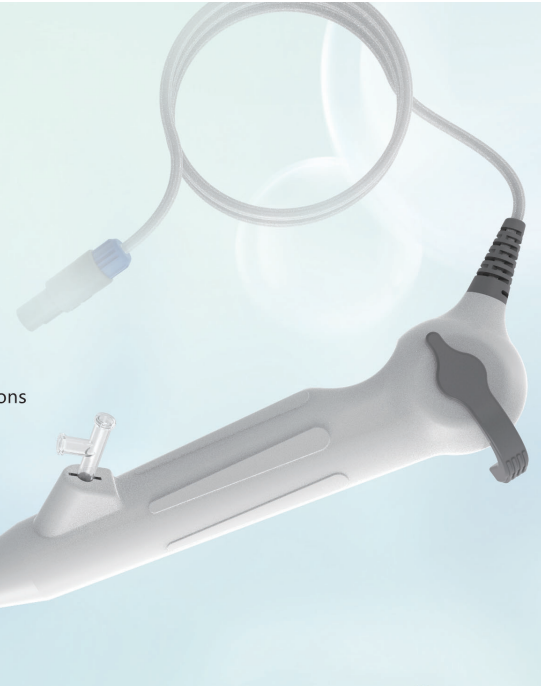

9.0Fr distal tip



- Excellent maneuverability
- Maximal deflection of 180° in two directions
- Self-locking deflection angle


- Bullet-like outline design for easy insertion

Opening of working channel

Lighting source


Lighting source

Video camera

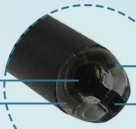

- Image resolution of 400×400 pixels
- Multilevel illumination adjustment


13.8Fr insertion tube outer diameter

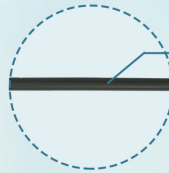

- The ultra-thin insertion tube facilitates operations and minimizes patient injury risks.
- 1:1 implementation of hand movements through to the distal end

Accessory port

Irrigation port

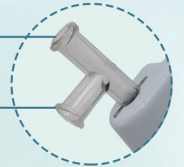

- Allows passage of accessory instruments and rapid fluid injection respectively during procedures.